


MEDICATION DISPENSER ACTIVATION GUIDE

CELLULAR AUTOMATED VERSION

USE YOUR UNIT ID AND THE PASSWORD: _____ TO LOG IN AT WWW.MEDREADY.NET.

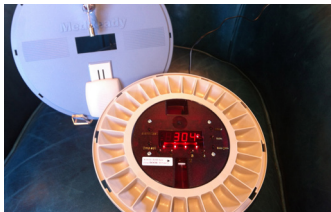

With the key, unlock the dispenser. Remove the lid by gently pressing down on the blue lid and turning it clockwise.




On the back of the dispenser, locate back-up battery compartment door. Open and insert the battery by plugging into small hole. To verify proper installation, the colon light will blink on the front panel. Plug the AC power supply into the socket on the dispenser AFTER connecting the battery. Plug the AC cord into an electrical outlet. Wait 20 minutes until the correct time of day appears on the clock face. If during this wait, an alarm sounds, slide the White T-BAR up until it clicks to silence the alarm.



Deactivate alarms 1,2,3, and 4 by pressing and holding the numbered buttons one at a time. The red light above each button will turn off.



Test the dispenser's ability to call in a medication dispenser alert. Press and hold the ALARM test button while you adjust the number on the clock from :30 to :03 minutes with the MIN button. This reduces the amount of time it takes to report an alert. Using alarm button 1, set ALARM 1 for a time 5 minutes past the displayed time on the clock. Once the alarm is set, make sure it is active by verifying the LED light is on above ALARM 1.

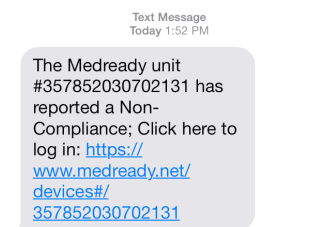

Wait for the alarm to activate. When it does, you will see the medication wheel turn and hear the alarm ring. Let it ring for the full three minutes. At the end of the time, you should see the large LED light come on. This means it is sending the notification signal. At this time you should receive the desired notification (automated phone call, text message or email).



Reset the notification time back to :30 or desired time by pressing and holding the Alarm Test button to display the current setting of :03 and using the MIN button to adjust the time back to :30 minutes (or desired time).



Set your dispense alarm times, using only the number of alarms you need. Leave the remaining alarms deactivated to avoid overdosing. To set the alarms, press and hold the desired alarm button (1-4). You should see the time change on the clock panel. Continue to hold the alarm set button as you change the time with the Hour, Min, and AM-PM buttons. The light must remain on for the alarm to be active.

ACTIVATION GUIDE CONT.



Fill the pill tray. **DO NOT FILL WITH NITROGLYCERINE**. Start with the slot that is to the right of the cup found directly below the white T-bar. This cup on the right should be filled with the NEXT upcoming dose. The remaining doses should follow around the wheel counter clockwise. The cup below the white T-bar should remain empty. Before putting the lid on the dispenser, make sure the white T-bar is slid all the way up in the locked position.



Place the lid on the medication dispenser. Make sure the door on the lid is closed as you place the lid on the base. Lock the lid with the key, remove the key, and return it to its secure location.

Benefits of a Cellular Medication Dispenser

- Manage Medications Online: Set dispensing times, change device settings, trigger early doses, view detailed reports, manage contacts.
- Notifications: Automated phone call, text message and/or email notifications available. Events such as the patient does not access their medication, the dispenser does not rotate correctly, or when the A/C power has been disconnected for 24 hours can trigger a notification.
- Tracking Your Remaining Doses: If you would like to track online the remaining doses available, when refilling the unit hold down the Time Set button and Alarm 1 at the same time until the clock face flashes. This resets the dose tracker on the website.

Navigating the Website

Once you are logged in, the changes you make on the website will be pushed to the device automatically (assuming the dispenser is set up and connected to power). The changes you make on the device (setting or activating alarms, etc.) will automatically be sent to the website. The website is simple to navigate and is presented in several sections:

- **Summary:** View the basic information and current status of the device.
- **Notifications:** View the current notifications. Email, text, and automated phone call notifications are available.
- **Alarms:** View current alarm settings. You can also send an alarm test, activate an early dose, or cancel an early dose.
- **Schedule:** See the actual compliance data and the remaining doses available.
- **Log:** Data tracking for all events that may occur (Alarm set, time set, loss of power, etc.)

The screenshot shows the MedReady website interface. At the top, there is a 'Go Back' link and a circled 'Edit Device and User' link with the text 'All times in unit timezone.' below it. A red arrow points to the 'Edit Device and User' link with the text 'Make changes here'. Below this, there is a list of settings: 'Off', 'Active', '3', 'd? Yes', 'd? Yes', and 'Yes'. On the right side, there is a sidebar with a menu containing 'Summary', 'Notifications', 'Alarms', 'Schedule', and 'Log'. Below the menu, there is a section for 'Cellular Unit ID: 013108001' and 'Time Zone: Eastern Time (US & C)'. Below that, there is an 'Information' section with the following details: Unit ID: 01310800161924, Model: Cellular, Organization: MedReady / Simp, Time Zone: Eastern Time (US, SIM: 89014103279388, and SIM Status: Suspended. A red arrow points to the 'Log' menu item with the text 'Use tabs to view current status'.

